Advertencia

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Warning

- Never under any circumstances attempt to open or inspect the internal components or accessories of the machine. If an inspection becomes necessary, please contact your supplier or a suitably competent technician.
- Never use the unit in humid or poorly ventilated environments, such as wet rooms, bathrooms etc. to prevent damage to the electrical components of your unit.
- Never use the unit outdoors; exposure to the elements can cause damage.
- If the power cord was destroyed, please ask a professional to service and replace it to avoid hazards, or contact the dealer for service.

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EC-Declaration of compliance

This product meets the following requirements:

1. Conditions set out in the Low Voltage Directive: 2014/35/EU.
2. Conditions set out in the Electromagnetic Compatibility Directive: 2014/30/EU.

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Recycling

Whenever any of our electrical products or their presentation packaging bear this symbol, it means that these products should not be treated as conventional household waste in Europe. To ensure proper treatment of this waste, dispose of it in accordance with local laws or as required for the disposal of electrical equipment. This initiative will help preserve save natural resources and improve standards of environmental protection in regard to the treatment of electrical waste.

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Welcome!

Thank you very much for purchase. Please read the instructions carefully before operating the instrument. You can also ask a professional for guidance to make sure you are using the instrument correctly. We will accept no responsibility for any accidents that occur due to incorrect operation.

The company reserves the right to amend the details of the instruments with no prior notice. If any mistakes are found in these instructions, please feel free to send us an e-mail informing us so that we can correct them.

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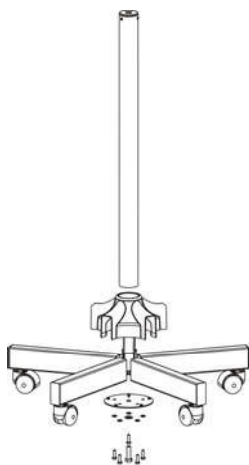
Instructions of use

1. Assemble the base according to the sketch shown below:
2. Insert the lamp into the base.
3. Turn on the switch.
4. Turn on the switch that is used to control the light intensity (the golden one) and adjust it to the proper intensity for the customer.
5. After the treatment, first adjust down the intensity, control switch to minimum, then switch off the red switch.

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Cautions

- Do not touch the bulb during the treatment or after within half an hour to avoid scalds.

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Maintenance

- Take good care of the bulb and clean it frequently. Do not use the lamp under unstable voltage.

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TECHNICAL SHEET

| | | |
|----------------------------|-------------------|-------------|
| Weights and sizes | Packing size (cm) | 77x30x24 |
| | Gross weight (kg) | 9 |
| | Net weight (kg) | 7,8 |
| Electrical characteristics | Rated voltage | 220V ~ 240V |
| | Frecuency | 50Hz/60Hz |
| | Rated power | 275 W |

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